


CARL ZEISS · JENA

Telegraphic Address: ZEISSWERKMED JENA

Branches: Berlin W. 9, Potsdamer Straße 139 / Hamburg, Rathausmarkt 8 / Vienna IX/3, Ferstelgasse 1
Buenos Aires, Bernardo de Irigoyen 330 / Tokyo, Tsukiji 33 a
Distributing Agencies: Milan, Via Lovanio 4 / Paris 5e, 26, rue Vauquelin / Madrid, Plaza de Canalejas 3, pr. d.
London W. C. 1, 8, Southampton Row / New York, 153 West 23rd Street

Prof. Koeppé's Eye Microscope

Like the Corneal Microscope, the Koeppé Eye Microscope is adapted for the stereo-microscopic examination and investigation of all parts and layers of the living eye, up to the retina, with natural (including redless and monochromatic) focal light, as furnished by Gullstrand's slit lamp. In addition, it lends itself to observation with *polarised focal light* and to important *ultra-microscopic* observations.

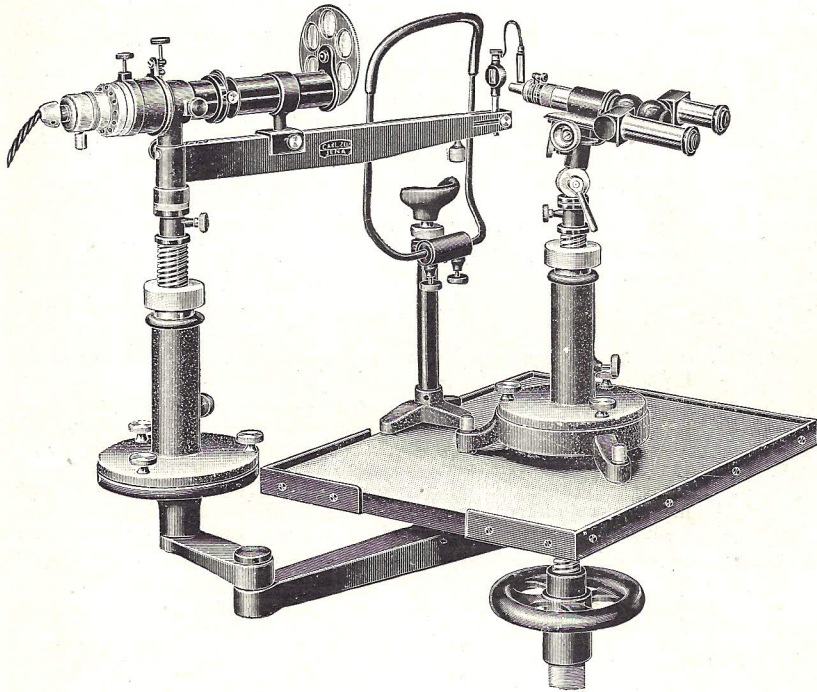

7532

The Eye Microscope (with polarising attachment) on stand and foot, together with the slit lamp (with screening tube, slit lamp lens and silvered mirror), all mounted on instrument table with jointed bracket arm, chin and head rests.

A single objective furnishes stereoscopic vision by a binocular microscope attachment. The latter may be arranged for observation by erect or inverted images. The latter arrangement is less costly, but the correct interpretation of

